


Cordis[™]

BRITE TIP[™]

Sheath Introducers

Illuminating the Path to
Precision Placement

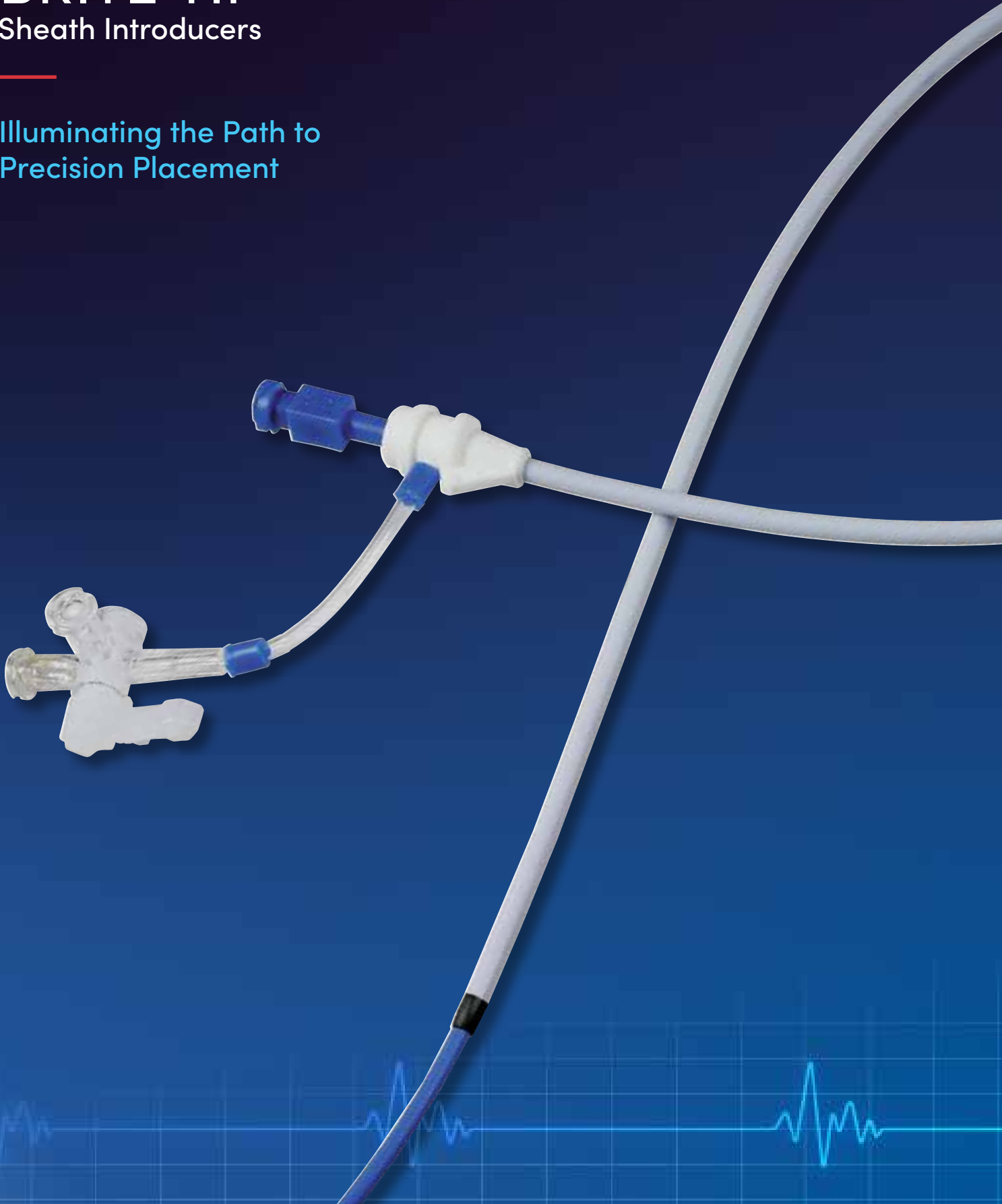

Kink resistant cannula

RESILIENT DUAL LAYER CANNULA


- Flexibility and support to access tortuous vessels

REINFORCED CANNULA WALL

- Kink resistance

SILICONE COATED CANNULA AND VESSEL DILATOR


- Lubricious, anti-thrombogenic vessel entry

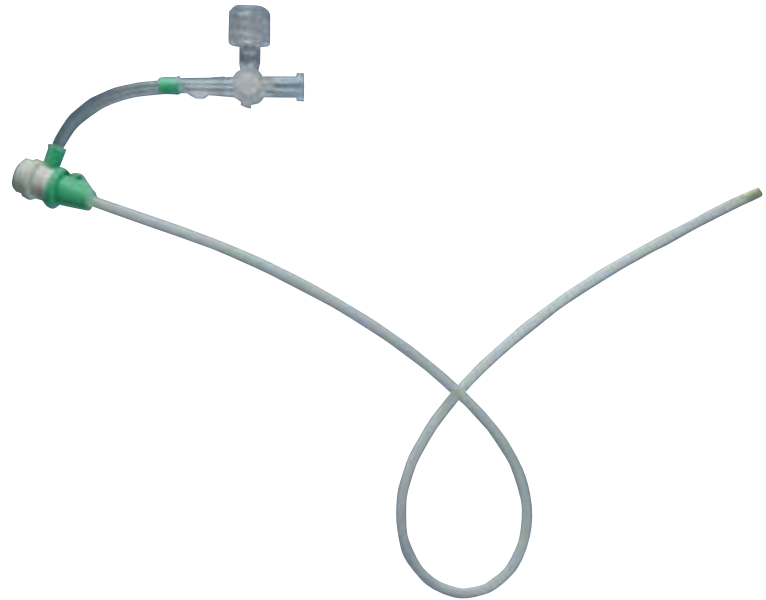

SLIX valve

SIX-CUT SPIRAL DESIGN FOR HEMOSTATIC VALVE

- Free catheter movement without blood reflux or air aspiration

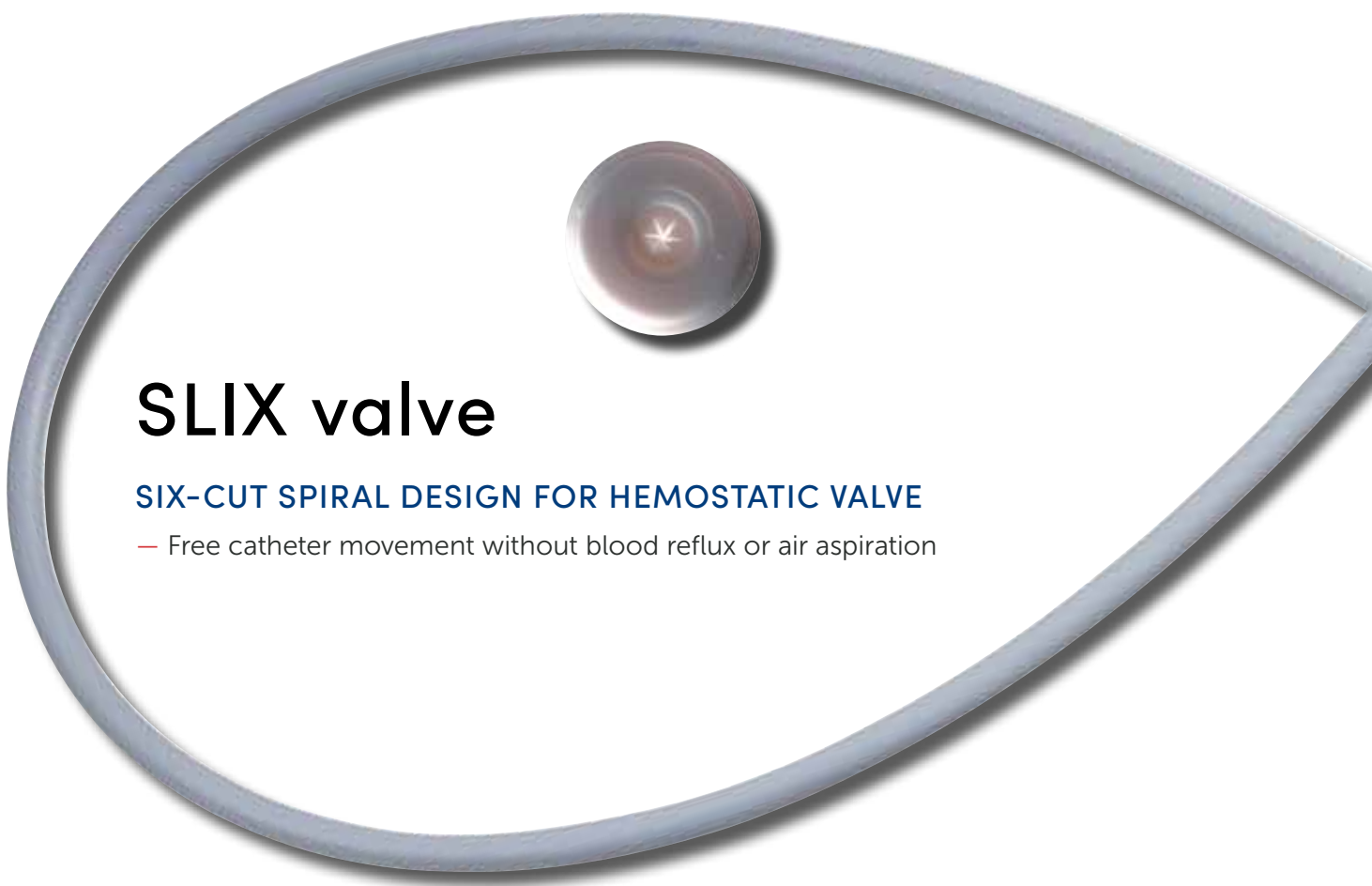


BRITE TIP™

Precise, atraumatic placement

Manicured tip and dilator design

TAPERED TRANSITIONS ON SHEATH AND DILATOR FOR:

- Atraumatic vessel entry
- Improved patient comfort

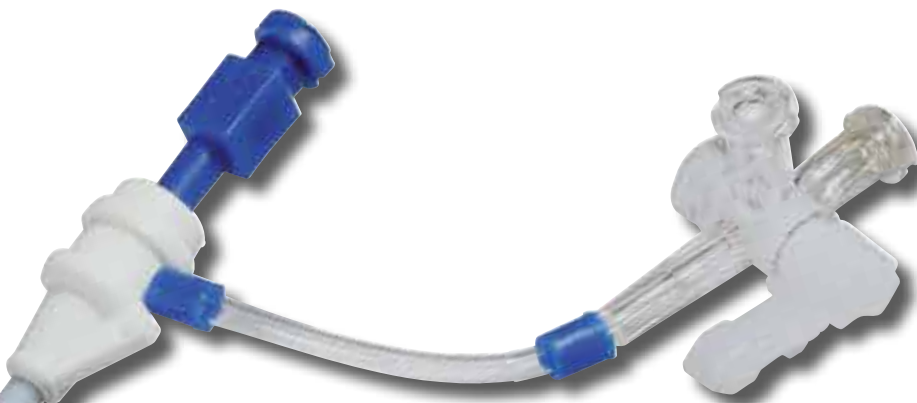

Radiopaque BRITE TIP™ Sheath

BLACK TUNGSTEN-FILLED DISTAL TIP

- Excellent visualization facilitates exact location of sheath's distal tip for accurate interventional positioning



BRITE TIP™

Sheath Introducers

ORDERING INFORMATION

Longer lengths available for a wide variety of interventional procedures



CATALOG NUMBERS

Cannula Length			4.5cm	5.5cm	11cm	23cm	35cm	45cm	55cm	70cm	90cm
French Size	Color Code	Dilator Tip (mm)	With Mini-Wire	With Mini-Wire	With Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire	Without Mini-Wire
4	■	25	SRA0116		401411M	401423M		SRAM17661			
5	■	25	SRA0117	SRA0035	401511M	401523M	401535M	401545M	401555M		
6	■	35	SRA0115	401605M	401611M	401623M	401635M	401645M	401655M	SRA0217	401690M
7	■	35		401705M	401711M	401723M	401735M	401745M	401755M		401790M
8	■	35		401805M	401811M	401823M	401835M	401845M	401855M		
9	■	35			401911M	401923M	401935M				
10	■	45			401011M	401023M					
11	■	45			401111M	401123M					

*Each catalog number must be ordered by the box

QUANTITY

- 5 units per package for 4.5cm, 5.5cm, 11cm and 23cm lengths
- 2 units per package for 35cm or longer lengths

For Healthcare Professional Only. Important information: Prior to use, refer to the instructions for use supplied with this device for indications, contraindications, side effects, suggested procedure, warnings and precautions. As part of its continuous product development policy, Cordis reserves the right to change product specifications without prior notification.

Please contact your Cordis representative for additional product availability information.

CORDIS, Cordis LOGO and BRITE TIP are trademarks of Cordis and may be registered in the US and/or in other countries. © 2022 Cordis. All Rights Reserved. 100546307-2 02/2022

